

Bath | water bath

digital, 6 / 11 / 22 l, up to 100°C ±0.1°C, flat lid included




MEGA LAB

Tel. 2310855844, 2106452848

Email. contact@megalab.gr

www.megalab.gr

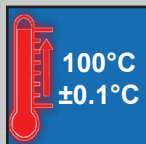

IQ | OQ
available
page 12 for details

WB 100°C, digital

Ideal for:

- water bath for general purpose

Features:

- user's self-compensation function to control difference between real temperature and indicated temperature (±10.0°C)
- stainless steel bath & lid for superior durability & high thermal efficiency / consistency, dome lid optional
- digital timer function (99hr 59min): delay time and operation time can be controlled
- storage function for the set values of temperature and timer
- backlight LCD display
- CE-certified and unique serial number for tracing

Safety mechanism:

- locking mode
- over-temperature and over-current protection
- sensor error protection

Controller:


- digital fuzzy controller with **Jog-Shuttle-switch** (turn + push)



WB-11 with flat lid (included)

Model	WB-6	WB-11	WB-22
Capacity	6 l	11 l	22 l
Temperature range & accuracy	room temperature +5°C - 100°C, ±0.1°C		
Temperature resolution	0.1°C-display, 0.1°C-control		
Temperature uniformity & probe	±1.0°C, PT100		
Heating power	800 W	1 kW	2 kW
Timer & alarm	99hr 59 min (delay & continuous run), error status & timer end		
Drain valve	-		mounted (Ø9.5 - 11 mm)
Effective space	300 x 155 x 110 mm	302 x 240 x 110 mm	500 x 295 x 110 mm
Internal dimensions (W x D x H)	300 x 155 x 150 mm	302 x 240 x 150 mm	500 x 295 x 150 mm
External dimensions (W x D x H)	325 x 180 x 260 mm	325 x 265 x 260 mm	535 x 330 x 260 mm
Packing size & gross weight	350 x 210 x 310 mm, 8 kg	350 x 300 x 310 mm, 10 kg	570 x 360 x 310 mm, 14 kg
Power supply	1 Phase AC 120V, 60 Hz or AC 230V, 50/60 Hz		
Order number 230V	DH.WB000106	DH.WB000111	DH.WB000122
Order number 120V	DH.WB000206	DH.WB000211	DH.WB000222

Accessories:

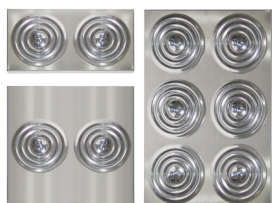

Lid with concentric ring, stainless steel

Measurements	Article	Order number
for WB-6, 2 rings, 330x182xh20 mm	BFC106	DH.WB030106
for WB-11, 2 rings, 329x269xh20 mm	BFC111	DH.WB030111
for WB-22, 6 rings, 539x334xh20 mm	BFC122	DH.WB030122



Dome lid for water bath WB, stainless steel

Measurements	Article	Order number
for WB-6, 324x176xh100 mm	BD206	DH.WB010206
for WB-11, 324x263xh100 mm	BD211	DH.WB010211
for WB-22, 533x328xh100 mm	BD222	DH.WB010222



Tube rack, stainless steel

Measurements	Article	Order number
for WB-6, 36x Ø17 mm, 220x125xh79 mm	BTH106	DH.WB020106
for WB-6, 15x Ø32 mm, 220x125xh79 mm	BTH206	DH.WB020206
for WB-11, 63x Ø17 mm, 220x205xh79 mm	BTH111	DH.WB020111
for WB-11, 25x Ø32 mm, 220x205xh79 mm	BTH211	DH.WB020211

for WB-22 please choose between 4x BTH106/206 or 2x BTH111/211